

This Page Is Inserted by IFW Operations  
and is not a part of the Official Record

## **BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning documents *will not* correct images,  
please do not report the images to the  
Image Problems Mailbox.**



## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

|  |                  |   |
|--|------------------|---|
| <b>(51) International Patent Classification:</b><br><b>G01N 33/74</b>  | <b>A1</b>        | <b>(11) International Publication Number:</b><br><b>WO 00/42437</b><br><b>(43) International Publication Date:</b><br>20 July 2000 (20.07.2000) |
| <b>(21) International Application Number:</b> PCT/US00/00855<br><b>(22) International Filing Date:</b> 13 January 2000 (13.01.2000)<br><b>(30) Priority Data:</b><br>09/231,422 14 January 1999 (14.01.1999) US<br>09/344,639 26 June 1999 (26.06.1999) US<br><b>(60) Parent Application or Grant</b><br>SCANTIBODIES LABORATORY, INC. [/]; O. VOYCE,<br>Brian, D.; O.   | <b>Published</b> |   |
| <b>(54) Title: METHODS FOR DIFFERENTIATING AND MONITORING PARATHYROID AND BONE STATUS RELATED DISEASES</b><br><b>(54) Titre: PROCEDE POUR DIFFERENTIER ET SURVEILLER LES MALADIES LIEES A L'ETAT DE LA PARATHYROIDE ET DES OS</b><br><br><b>(57) Abstract</b><br><p>The present invention relates to novel methods and devices for differentiating in a patient parathyroid diseases, such as hyperparathyroidism and related bone diseases, from normal or non-disease states. One detects whole or non-fragmented (1 to 84) parathyroid hormone in a biological sample and also a large non-whole parathyroid hormone peptide fragment that can function as a parathyroid hormone antagonist. By either comparing values or using independently the value of either the large non-whole parathyroid hormone peptide fragment, the whole parathyroid hormone, or the combination of these values one is able to differentiate parathyroid and bone related disease states, as well as differentiate such states from normal states.</p> <b>(57) Abrégé</b><br><p>L'invention concerne de nouveaux procédés et dispositifs pour distinguer chez un sujet des maladies de la parathyroïde, telles que l'hyperparathyroïdie et les maladies osseuses associées, des états normaux ou de l'absence de maladie. On détecte dans un prélèvement biologique l'hormone parathyroïde entière ou non fragmentée (1 à 84) ainsi qu'un grand fragment non entier du peptide de l'hormone parathyroïde qui peut fonctionner comme un antagoniste de l'hormone parathyroïde. En comparant les valeurs ou en utilisant indépendamment la valeur du grand fragment non entier du peptide de l'hormone parathyroïde, de l'hormone parathyroïde entière et de la combinaison de ces valeurs, on peut faire la distinction entre les états pathologiques liés à la parathyroïde ou aux os de même qu'entre ces états pathologiques et les états normaux.</p> |                  |   |

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

|  |           |   |
|--|-----------|---|
| <b>(51) International Patent Classification 7 :</b><br><b>G01N 33/74</b>   | <b>A1</b> | <b>(11) International Publication Number:</b> <b>WO 00/42437</b><br><b>(43) International Publication Date:</b> 20 July 2000 (20.07.00)   |
| <b>(21) International Application Number:</b> PCT/US00/00855<br><b>(22) International Filing Date:</b> 13 January 2000 (13.01.00)<br><br><b>(30) Priority Data:</b><br>09/231,422 14 January 1999 (14.01.99) US<br>09/344,639 26 June 1999 (26.06.99) US<br><br><b>(71) Applicant:</b> SCANTIBODIES LABORATORY, INC. [US/US];<br>9336 Abraham Way, Santee, CA 92071 (US).<br><br><b>(74) Agent:</b> VOYCE, Brian, D.; Suite C204, 1100 Possum Trot<br>Road, North Myrtle Beach, SC 29582 (US).   |           | <b>(81) Designated States:</b> CA, JP, European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).<br><br><b>Published</b><br><i>With international search report.</i><br><i>Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i> |
| <b>(54) Title:</b> METHODS FOR DIFFERENTIATING AND MONITORING PARATHYROID AND BONE STATUS RELATED DISEASES<br><br><b>(57) Abstract</b><br><br>The present invention relates to novel methods and devices for differentiating in a patient parathyroid diseases, such as hyperparathyroidism and related bone diseases, from normal or non-disease states. One detects whole or non-fragmented (1 to 84) parathyroid hormone in a biological sample and also a large non-whole parathyroid hormone peptide fragment that can function as a parathyroid hormone antagonist. By either comparing values or using independently the value of either the large non-whole parathyroid hormone peptide fragment, the whole parathyroid hormone, or the combination of these values one is able to differentiate parathyroid and bone related disease states, as well as differentiate such states from normal states. |           |   |

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

|    |                          |    |  |    |  |    |                          |
|----|--------------------------|----|--|----|--|----|--------------------------|
| AL | Albania                  | ES | Spain                                    | LS | Lesotho                                      | SI | Slovenia                 |
| AM | Armenia                  | FI | Finland                                  | LT | Lithuania                                    | SK | Slovakia                 |
| AT | Austria                  | FR | France                                   | LU | Luxembourg                                   | SN | Senegal                  |
| AU | Australia                | GA | Gabon                                    | LV | Latvia                                       | SZ | Swaziland                |
| AZ | Azerbaijan               | GB | United Kingdom                           | MC | Monaco                                       | TD | Chad                     |
| BA | Bosnia and Herzegovina   | GE | Georgia                                  | MD | Republic of Moldova                          | TG | Togo                     |
| BB | Barbados                 | GH | Ghana                                    | MG | Madagascar                                   | TJ | Tajikistan               |
| BE | Belgium                  | GN | Guinea                                   | MK | The former Yugoslav<br>Republic of Macedonia | TM | Turkmenistan             |
| BF | Burkina Faso             | GR | Greece                                   | ML | Mali   | TR | Turkey                   |
| BG | Bulgaria                 | HU | Hungary                                  | MN | Mongolia                                     | TT | Trinidad and Tobago      |
| BJ | Benin                    | IE | Ireland                                  | MR | Mauritania                                   | UA | Ukraine                  |
| BR | Brazil                   | IL | Israel                                   | MW | Malawi                                       | UG | Uganda                   |
| BY | Belarus                  | IS | Iceland                                  | MX | Mexico                                       | US | United States of America |
| CA | Canada                   | IT | Italy                                    | NE | Niger  | UZ | Uzbekistan               |
| CF | Central African Republic | JP | Japan                                    | NL | Netherlands                                  | VN | Viet Nam                 |
| CG | Congo                    | KE | Kenya                                    | NO | Norway                                       | YU | Yugoslavia               |
| CH | Switzerland              | KG | Kyrgyzstan                               | NZ | New Zealand                                  | ZW | Zimbabwe                 |
| CI | Côte d'Ivoire            | KP | Democratic People's<br>Republic of Korea | PL | Poland                                       |    |                          |
| CM | Cameroon                 | KR | Republic of Korea                        | PT | Portugal                                     |    |                          |
| CN | China                    | KZ | Kazakhstan                               | RO | Romania                                      |    |                          |
| CU | Cuba                     | LC | Saint Lucia                              | RU | Russian Federation                           |    |                          |
| CZ | Czech Republic           | LI | Liechtenstein                            | SD | Sudan  |    |                          |
| DE | Germany                  | LK | Sri Lanka                                | SE | Sweden                                       |    |                          |
| DK | Denmark                  | LR | Liberia                                  | SG | Singapore                                    |    |                          |
| EE | Estonia                  |    |  |    |  |    |                          |

**Description**

5

10

15

20

25

30

35

40

45

50

55

**METHODS FOR DIFFERENTIATING AND MONITORING  
PARATHYROID AND BONE STATUS RELATED DISEASES**

**TECHNICAL FIELD**

The present invention relates to novel methods and devices for differentiating in a patient parathyroid diseases, such as hyperparathyroidism, from normal or non-disease states. One detects whole or non-fragmented (1 to 84) parathyroid hormone in a biological sample and also a large non-whole parathyroid hormone peptide fragment that can function as a parathyroid hormone antagonist. By either comparing values or using independently the value of either the large non-whole parathyroid hormone peptide fragment, the whole parathyroid hormone, or the combination of these values one can differentiate parathyroid and bone related disease states, as well as differentiate such states from normal states.

**RELATED APPLICATIONS**

The present application is a continuation-in-part of a non-provisional utility patent application filed in the United States Patent and Trademark Office, Serial Number 08/231,422.

**BACKGROUND ART**

Calcium plays an indispensable role in cell permeability, the formation of bones and teeth, blood coagulation, transmission of nerve impulse, and normal muscle contraction. The concentration of calcium ions in the blood is, along with calcitriol and calcitonin, regulated mainly by parathyroid hormone (PTH). Although calcium intake and excretion may vary, PTH serves through a feedback mechanism to maintain a steady concentration

5 of calcium in cells and surrounding fluids. When serum calcium lowers, the parathyroid glands secrete PTH, affecting the release of stored calcium. When serum calcium increases, stored calcium release is retarded through lowered secretions of PTH.

10 5 The complete form of human PTH, sometimes referred to in the art as hPTH but referred to in the present invention either as whole PTH or wPTH, is a unique 84 amino acid peptide (SEQ ID NO. 1), as is shown in FIGURE 1. Researchers have found that this  
15 peptide has an anabolic effect on bone that involves a domain for protein kinase C activation (amino acid residues 28 to 34) as well as a domain for adenylate cyclase  
20 activation (amino acid residues 1 to 7). However, various catabolic forms of clipped or fragmented PTH peptides also are found in circulation, most likely formed by intraglandular or peripheral metabolism. For example, whole PTH can be cleaved between amino acids 34 and 35 to produce a (1-34) PTH N-terminal fragment and a (35-84) PTH C-terminal fragment. Likewise, clipping can occur between either amino acids 36 and 37  
25 15 or 37 and 38. Recently, a large PTH fragment referred to as "non-(1-84) PTH" has been disclosed which is clipped closer to the N-terminal end of PTH. (See R. LePage *et alia*, "A non-(1-84) circulating parathyroid hormone (PTH) fragment interferes significantly with intact PTH commercial assay measurements in uremic samples" Clin Chem (1998); 44: 805-810.)

20 20 The clinical need for accurate measurement of PTH is well demonstrated. Serum  
35 PTH level is one of the most important indices for patients with the following diseases: familial hypocalciuria; hypercalcemia; multiple endocrine neoplasia types I and II; osteoporosis, Paget's bone disease; primary hyperparathyroidism - caused by primary  
40 25 hyperplasia or adenoma of the parathyroid glands; pseudohypoparathyroidism; and renal failure, which can cause secondary hyperparathyroidism.

45 PTH plays a role in the course of disease in a patient with chronic renal failure. Renal osteodystrophy (RO) is a complex skeletal disease comprising osteitis fibrosa  
30 cystica (caused by PTH excess), osteomalacia - unmineralized bone matrix (caused by

5 vitamin D deficiency), extraskeletal calcification/ossification (caused by abnormal calcium  
and phosphorus metabolism), and adynamic bone disease (contributed to by PTH  
suppression). Chronic renal failure patients can develop RO. Failing kidneys increase  
10 serum phosphorus (hyperphosphoremia) and decrease 1,25-dihydroxyvitamin D (1,25-D)  
5 production by the kidney. The former results in secondary hyperparathyroidism from  
decreased gastrointestinal calcium absorption and osteitis fibrosa cystica from increased  
PTH in response to an increase in serum phosphorus. The later causes hypocalcemia and  
15 osteomalacia. With the onset of secondary hyperparathyroidism, the parathyroid gland  
becomes less responsive to its hormonal regulators because of decreased expression of its  
10 calcium and vitamin D receptors. Serum calcium drops. RO can lead to digital gangrene,  
bone pain, bone fractures, and muscle weakness.

Determining circulating biologically active PTH levels in humans has been  
challenging. One major problem is that PTH is found at low levels, normally 10pg/mL to  
25 65 pg/mL. Coupled with extremely low circulating levels is the problem of the  
15 heterogeneity of PTH and its many circulating fragments. In many cases, immunoassays  
have faced substantial and significant interference from circulating PTH fragments. For  
example, some commercially available PTH kits have almost 100% cross-reactivity with  
30 the non-(1-84) PTH fragment, (see the LePage article).

20 PTH immunoassays have varied over the years. One early approach is a double  
antibody precipitation immunoassay found in U. S. 4,369,138 to Arnold W. Lindall *et alia*.  
35 A first antibody has a high affinity for a (65-84) PTH fragment. A radioactive labeled (65-  
84) PTH peptide is added to the sample with the first antibody to compete for the  
40 25 endogenous unlabeled peptide. A second antibody is added which binds to any first  
antibody and radioactive labeled PTH fragment complex, thereby forming a precipitate.  
Both precipitate and supernatant can be measured for radioactive activity, and endogenous  
45 PTH levels can be calculated therefrom.



5 In an effort to overcome PTH fragment interference, immunoradiometric two-site  
assays for intact PTH (I-PTH) have been introduced, such as Allegro® Intact PTH assay  
10 by the Nichol's Institute of San Juan Capistrano, California. In one version, a capture  
5 antibody specifically binds to the C-terminal portion of hPTH while a labeled antibody  
specifically binds to the N-terminal portion of the captured hPTH. In another, two  
monoclonal antibodies were used, both of which attached to the N-terminal portion of  
15 hPTH. Unfortunately, these assays have problems in that they measure but do not  
discriminate between wPTH and non-whole PTH peptide fragments. This inability comes  
10 to the fore in hyperparathyroid patients and renal failure patients who have significant  
endogenous concentrations of large, non-whole PTH fragments.

20 Recently, researchers have made a specific binding assay directed to the large N-  
terminal PTH fragments. (See. Gao, Ping *et alia* "Immunochemiluminometric assay  
25 with two monoclonal antibodies against the N-terminal sequence of human parathyroid  
hormone", Clinica Chimica Acta 245 (1996) 39-59.) This immunochemiluminometric  
assay uses two monoclonal antibodies to detect N-terminal (1-34) PTH fragments but not  
30 mid-portion PTH fragments or C-terminal PTH fragments. A key factor in the design of  
these assays is to eliminate any reaction with C-terminal PTH fragments.

#### 35 DISCLOSURE OF THE INVENTION

40 The present invention relates to novel methods and devices for differentiating in a  
25 patient parathyroid diseases, (such as primary hyperparathyroidism, secondary  
hyperparathyroidism, and stages thereof), from normal or non-disease states; for  
monitoring the function of parathyroid glands either during or after treatment, i.e., intra-  
45 operation and after operation parathyroid function monitoring as well as therapeutic  
treatment; and also for monitoring the effects of therapeutic treatments for parathyroid

5 related bone diseases and hyperparathyroidism. One detects the level in the serum or  
blood of at least one of three different parameters, namely, whole or non-fragmented  
parathyroid hormone in a biological sample, a large non-whole parathyroid hormone  
10 peptide fragment that can function as a parathyroid hormone antagonist, or the  
5 combination of the two values. By comparing the two values or by examining  
independently one of the above three values, one can differentiate parathyroid and bone  
disease states, as well as differentiate such states from normal states, as the relationship  
15 between these values, as well as the values themselves, change significantly between a  
normal person and a patient with a parathyroid disease.

10  
20 The present invention incorporates a discovery that a large, non-whole  
PTH peptide fragment, a peptide having an amino acid sequence from between (SEQ ID  
No.2 [PTH<sub>1-34</sub>]) and (SEQ ID No. 3 [PTH<sub>34-14</sub>]), functions *in vivo* as a wPTH antagonist  
or inhibitor (PIN), (see FIGURE 12). In other words, the binding of wPTH to PTH  
25 15 receptors and the subsequent biological activity are affected by the presence of this PIN  
peptide fragment. The PTH receptors can be tied up with respect to PTH or PTH  
analogs in that the PTH binding site is blocked. The relationship between the  
30 concentrations of wPTH and PIN vary with PTH related disease states, and thus, are  
indicative of such states. Equally useful in view of the discovery of the antagonist nature  
20 of PIN, the present invention relates to novel methods and devices for monitoring  
parathyroid related bone diseases, and resultant bone loss or build-up. Increased  
35 amounts of PIN can inhibit the calcium releasing activity of PTH.

40 25 In making a measurement of wPTH, one does not want to detect PIN. The  
method for measuring the amount of wPTH in a sample such as serum, plasma, or blood  
comprises four general steps which can vary depending upon whether one uses a first  
antibody or antibody fragment specific for the PTH peptide SER-VAL-SER-GLU-ILE-  
45 GLN-LEU-MET (SEQ ID No. 4), wherein at least four amino acids are part of the  
antibody reactive portion of the peptide either as a signal antibody or a capture antibody in

5 conventional immunoassay formats. (One can also use an analogous peptide present in  
other species, such as a rat peptide in which the first amino acid serine is substituted with  
an alanine.) Used either as a signal antibody or as a capture antibody, enough antibody is  
10 added to bind all wPTH present. Next, one allows the first antibody to bind to any wPTH  
5 present, thereby forming a complex. A specific binding label comprised of a second  
antibody and a conventional immunoassay label, such as chemiluminescent agents,  
colorimetric agents, energy transfer agents, enzymes, fluorescent agents, and  
15 radioisotopes, is used to label the complex, preferably at the C-terminal end of wPTH, and  
can be added either substantially simultaneously with the first antibody or subsequent  
10 thereto. Finally, one uses conventional techniques to measure the amount of labeled  
complex, and thereby calculate wPTH levels in the sample. If used as a signal antibody,  
20 then the first antibody still attaches at the N-terminal end, but the second antibody would  
serve as a capture antibody that attaches at the C-terminal end.

25 15 In making a measurement of PIN, one can either measure it directly, or indirectly.  
An indirect measurement can be made by first measuring wPTH and then measuring total  
PTH. Subtracting the wPTH value from the total PTH value, one derives the PIN value.  
30 (For the purposes of the present invention, "total PTH" refers to the sum of wPTH, the  
naturally occurring predominant PTH receptor binding agonist, and PIN, the naturally  
20 occurring predominant PTH receptor binding antagonist.) A total PTH assay detects both  
PIN and wPTH by detecting the N-terminal end of PTH not at SEQ ID No. 4, the very  
35 end of the N-terminal. By detecting between about amino acids 7 to 38 of PTH, the assay  
can detect both. A commercially available assay for total PTH is available from  
Scantibodies Laboratory, Inc. of Santee, California. A direct measurement of total PTH  
40 25 can be made by using an antibody or antibody fragment specific for a portion of the PTH  
peptide LEU-MET-HIS-ASN-LEU-GLY-LYS-HIS-LEU-ALA-SER-VAL-GLU-ARG-  
MET-GLN-TRP-LEU-ARG-LYS-LYS-LEU-GLN-ASP-VAL-HIS-ASN-PHE-VAL-  
45 ALA-LEU-GLY (SEQ ID No. 5), which comprises amino acids 7 to 38 of PTH,  
(preferably between amino acids 9 to 34), wherein at least four amino acids are part of the

antibody reactive portion of the peptide. Such an antibody or antibody fragment can be used in conventional immunoassay formats either as a signal antibody or a capture antibody.

To differentiate between parathyroid disease states and the normal state or to monitor the effects of therapeutic treatment for parathyroid disease states, one can compare the relationship between the values of wPTH, PIN, or total PTH, (the combination of wPTH and PIN), in other words, the relationship between the values of PIN and total PTH, between PIN and whole PTH, or between whole PTH and total PTH. For example, one can use a proportion between wPTH and total PTH, between PIN and total PTH, or between PIN and wPTH. (Comparisons can even take the form of a neural network of all these factors.) Regardless of the comparative method chosen, these values change significantly between a normal person and a patient with a parathyroid disease and between various stages of parathyroid diseases.

Alternatively, one can either differentiate between parathyroid disease states and the normal state or monitor the effects of therapeutic treatment for parathyroid disease states by examining independently the value of either wPTH, PIN, or total PTH alone.

#### **BRIEF DESCRIPTION OF THE DRAWINGS**

FIGURE 1 is a diagrammatic view of human wPTH.

FIGURE 2 is a diagrammatic view of a wPTH assay using the present antibody as a tracer element.

FIGURE 3 is a diagrammatic view of a wPTH assay using the present antibody as a capture element.

5  
FIGURE 4 is a graph showing a standard curve for a wPTH assay.

10  
5  
FIGURE 5 is a graph comparing a conventional I-PTH assay with the present wPTH assay for healthy normal persons with "normal" PTH values.

15  
FIGURE 6 is a diagrammatic view showing interference from PIN in conventional I-PTH assays.

10  
20  
FIGURE 7 is a graph comparing a conventional I-PTH assay with the present wPTH assay for patients with chronic uremia.

25  
15  
FIGURE 8 is a graph showing the distribution of wPTH values for healthy normal persons, patients with primary hyperparathyroidism, and patients with chronic uremia.

30  
20  
FIGURE 9 is a diagrammatic view showing how PIN blocks the action of wPTH at the receptor level, thereby making the person insensitive to the biological effects of wPTH.

35  
FIGURE 10 is a graph demonstrating complete cross-reactivity of wPTH and PIN in a total PTH assay used in the present invention.

40  
25  
FIGURE 11 is a graph demonstrating how the whole PTH assay used in the present invention does not detect to PIN.

45  
FIGURE 12 is a graph demonstrating how PIN is an *in vivo* inhibitor of wPTH.

5

**BEST MODES FOR CARRYING OUT THE INVENTION**

10

5 In disclosing the present invention, one should remember that there are a number of closely analogous, species dependent forms of PTH. The amino acid sequence of hPTH is shown in FIGURE 1. However, for rat PTH, bovine PTH, or porcine PTH, for example, one finds the substitutions at some of the amino acids in the hPTH sequence. For the purposes of the present invention, one can use interchangeably antibodies or antibody fragments to forms of these PTHs, although it is preferred to use an antibody with specificity for PTH having a sequence matching the species in which the PTH measurements are made.

20

**Whole PTH immunoassay**

25

15 A preferred embodiment of the present invention is an immunoradiometric assay (IRMA), often referred to as a sandwich assay, as shown FIGURES 2 and 3. Elements employed in such an assay (10) include a capture antibody (12) attached to a solid support (14) and a signal antibody (16) having a label (18), attached thereto (20). Typically, one selects a capture antibody that is specific for C-terminal PTH fragments (22), while the label antibody is specific for the initial wPTH peptide sequence which comprises a domain for adenylate cyclase activation (24), as shown in FIGURE 2. However, one could reverse the specificity of these antibodies, as is shown in FIGURE 3.

35

40

25 Alternatively, one could create an immunoassay in which wPTH is either precipitated from solution or otherwise differentiated in a solution, as in conventional precipitating assays or turbidometric assays. For example, one can use at least three antibodies to form a precipitating mass. In addition to the initial wPTH sequence antibody and a C-terminal antibody, one can use at least a third antibody which attaches to the mid portion of PTH. The combined mass of wPTH and the at least three antibodies would form a labeled precipitating mass which can be measured by conventional techniques.

45

50

55

5 Another method would be to couple the initial wPTH sequence antibody to colloidal solid supports, such as latex particles.

10 More specifically, one can create a signal antibody by iodinating 50 micrograms of  
5 affinity purified goat anti-(1-6) PTH antibody (Scantibodies Laboratory, Inc., Santee California, U.S.A.) by oxidation with chloramine T, incubation for 25 seconds at room temperature with 1 millicurie of 125-I radioisotope and reduction with sodium  
15 metabisulfate. Unincorporated 125-I radioisotope is separated from the 125-I-Goat anti-(1-6) PTH signal antibody by, passing the iodination mixture over a PD-10 desalting  
10 column (Pharmacia, Uppsala, Sweden) and following the manufacturers instructions. The fractions collected from the desalting column are measured in a gamma counter and those  
20 fractions representing the 125-I-goat anti-(1-6) PTH antibody are pooled and diluted to approximately 300,000 DPM (disintegrations per minute) per 100 microliters. This solution is the tracer solution to be used in the whole PTH IRMA.

25 15 Capture antibody coated tubes can be created by attaching affinity purified goat anti PTH 39-84 antibody, (Scantibodies Laboratory, Inc., Santee, California, U.S.A.), to  
30 12 x 75 mm polystyrene tubes (Nunc, Denmark) by means of passive absorption techniques which are known to those of skill in the art. The tubes are emptied and dried,  
20 creating solid phase antibody coated tubes.

35 To conduct a whole PTH assay of a sample, 200 microliter samples of human serum are added to the solid phase antibody coated tubes. To each tube is added 100  
40 25 microliters of the tracer solution (labeled goat anti-(1-6) PTH signal antibody). The tubes are incubated at room temperature with shaking at 170 rpm for 20-22 hours. During this time the immunochemical reaction of forming the sandwich of {solid phase goat anti-(39-84) PTH antibody} -- {whole PTH} -- {125-I-goat anti-(1-6) PTH antibody} takes place.  
45 Following this incubation, the test tubes are washed with distilled water. Radioactivity on the solid phase, which amount corresponds to the quantity of wPTH present, is measured

5 using a gamma counter. The radioactivity data for the samples is processed by  
conventional analysis with use of the results from standards and controls and a computer  
software in order that the concentration of whole PTH in the samples may be ascertained.  
10 FIGURE 4 shows a standard curve for such an assay.

5  
Initial whole PTH sequence peptide

15 In order to make the signal antibody in the above assay, first one makes a synthetic  
PTH peptide corresponding either to hPTH (Ser - Val - Ser - Glu - Ile - Gln - Leu - Met),  
rat PTH (Ala - Val - Ser - Glu - Ile - Gln - Leu - Met), or at least four amino acids in the  
10 common sequence. The selected peptide can play two roles in making an assay, first as a  
specific source for creating a polyclonal antibody or monoclonal antibody source for signal  
20 antibody or capture antibody, and second as part of an affinity purification means for  
isolating the desired signal antibody or capture antibody.

25 15 Briefly, such a peptide can be synthesized on an Applied Biosystems, Inc. (Foster  
City, California, U.S.A.) Model 431 automated peptide synthesizer employing Fmoc (9-  
fluoronylmethoxycarbonyl) as the alpha-amino protecting group. All amino acids and  
30 solvents are from Applied Biosystems and are of synthesis grade. Following synthesis, the  
peptide is cleaved from the resin, and side chains are de-blocked, using a cleavage cocktail  
20 containing 6.67% phenol, 4.4% (v/v) thioanisole and 8.8% ethanedithiol in trifluoroacetic  
acid (TFA). The cleaved peptide is precipitated and washed several times in cold diethyl  
35 ether. It is then dissolved in water and lyophilized. The crude peptide is subjected to  
amino acid analysis (Waters PICO-TAG System, Boston, Massachusetts, U.S.A.) and  
reversed-phase HPLC using a VYDAC (TM) C8 column with 0.1% TFA in water and  
40 25 99.9% acetonitrile in 0.1% TFA as the mobile buffers. The presence of a single major  
peak along with the appropriate amino acid composition is taken as evidence that the  
peptide is suitable for further use.



5           The resulting peptide is then attached to cross linked agarose beads (activated  
Sephacrose 4B from Pharmacia, Uppsala, Sweden) according to instructions from the  
manufacturer. Armed with the initial peptide sequence on a bead, one can affinity purify a  
10 polyclonal antibody serum source to isolate the initial sequence antibody for the wPTH  
5 immunoassay.

Initial sequence whole PTH antibody

15           To create an affinity-purified anti-(1-6) PTH antibody, one first uses a selected  
initial PTH sequence peptide as described above as part of an immunogen for injection  
10 into a goat. The peptide can be used either by itself as an injectible immunogen,  
incorporated into a non PTH peptide having a molecular weight, typically, of between  
20 about 5,000 and 10,000,000, or as part of the wPTH complete sequence. The  
immunogen is mixed with an equal volume of Freund's complete adjuvant which is a  
mixture of light mineral oil, Arlacel detergent, and inactivated mycobacterium tuberculosis  
25 15 bacilli. The resulting mixture is homogenized to produce an aqueous/oil emulsion which is  
injected into the animal (typically a goat) for the primary immunization. The immunogen  
dose is approximately 50-400 micrograms. The goats are injected monthly with the same  
dose of immunogen complex except no mycobacterium tuberculosis bacilli is used in these  
30 subsequent injections. The goats are bled monthly, approximately three months after the  
20 primary immunization. The serum (or antiserum) is derived from each bleeding by  
separating the red blood cells from the blood by centrifugation and removing the  
35 antiserum which is rich in (1-6) PTH antibodies.

40           To purify the antiserum for the desired (1-6) PTH antibody, one packs a separation  
25 column with the initial PTH sequence peptide bound beads described above, washes the  
column and equilibrates it with 0.01 M phosphate buffered saline (PBS). The antiserum  
is loaded onto the column and washed with 0.01 M PBS in order to remove antibodies  
45 without the (1-6) PTH specificity. The bound specific goat anti-(1-6) PTH polyclonal  
antibody is eluted from the solid phase PTH 1-6 in the column by passing an elution

5 solution of 0.1 M glycine hydrochloride buffer, pH 2.5 through the column. The eluted  
polyclonal antibody is neutralized after it leaves the column with either the addition of  
1.0 M phosphate buffer, pH 7.5 or by a buffer exchange with 0.01 M PBS, as is known  
10 to those of skill in the art. The polyclonal antibody is stored at 2-8 degrees centigrade.

Comparison between whole PTH and total PTH assays

15 The present wPTH IRMA assay was compared to a conventional intact PTH or I-  
10 PTH immunoassay, the Allegro Nichols Intact-PTH assay, (which is commercially  
available and made by Nichols Institute Diagnostics of San Juan Capistrano, California,  
U.S.A.), in both PTH normal persons and those suffering from chronic uremia. This I-  
20 PTH immunoassay, due to its 100% cross reactivity between PIN and wPTH, is in  
actuality a total PTH assay, (see FIGURE 10).

15  
25 FIGURE 5 shows the results for 34 normal human serum samples from healthy  
subjects which were assayed both by the present wPTH IRMA and the above I-PTH  
assay. In every case, the level of wPTH detected by the IRMA is lower than that reported  
by the I-PTH assay, demonstrating the ability of the present IRMA to avoid detecting the  
30 20 interfering large, non (1-84) PTH fragment detected by the I-PTH assay, (see FIGURE  
11). FIGURE 6 illustrates how such interference can occur. An N-terminal PTH specific  
signal antibody which is not specific to the initial PTH peptide sequence, as in the present  
35 invention, can detect not only wPTH (as in the upper part of FIGURE 6), but also can  
detect PIN, the large, non (1-84) PTH fragment, (as in the lower part of FIGURE 6).

25  
40 A comparison of assay results for 157 chronic uremic patients is shown in  
FIGURE 7. Serum samples from these patients were measured using the wPTH IRMA  
and the above I-PTH assay. In every case the wPTH levels are lower than I-PTH values.

Clinical Use

The present wPTH and PIN assays have been used in a clinical setting involving 188 persons. The group included 31 persons having normal healthy parathyroid glands and 157 patients with chronic uremia who are undergoing dialysis on a continuous basis. Each person had a blood sample drawn which was assayed using a wPTH assay from Scantibodies Laboratory, Inc. as well as an I-PTH assay from Nichols Institute which gave total PTH values.

Table 1 shows the results individually and comparatively, of the wPTH, PIN, and total PTH assays from chronic uremic patients on dialysis.

TABLE 1

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN<br>to<br>Total PTH | PIN<br>to<br>Whole PTH | Whole PTH<br>to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|------------------------|------------------------|------------------------------|
| 1           | 1410               | 740                | 670          | 48%                    | 91%                    | 52%                          |
| 2           | 185                | 89                 | 96           | 52%                    | 108%                   | 48%                          |
| 3           | 231                | 104                | 127          | 55%                    | 122%                   | 45%                          |
| 4           | 1020               | 590                | 430          | 42%                    | 73%                    | 53%                          |
| 5           | 270                | 159                | 111          | 41%                    | 70%                    | 59%                          |
| 6           | 201                | 100                | 101          | 50%                    | 101%                   | 50%                          |
| 7           | 380                | 100                | 280          | 74%                    | 280%                   | 26%                          |
| 8           | 460                | 277                | 183          | 40%                    | 66%                    | 60%                          |
| 9           | 380                | 197                | 183          | 48%                    | 93%                    | 52%                          |
| 10          | 880                | 522                | 358          | 41%                    | 69%                    | 59%                          |
| 11          | 310                | 154                | 156          | 50%                    | 101%                   | 50%                          |
| 12          | 880                | 451                | 429          | 49%                    | 95%                    | 51%                          |
| 13          | 670                | 418                | 252          | 38%                    | 60%                    | 63%                          |
| 14          | 390                | 221                | 169          | 43%                    | 76%                    | 57%                          |
| 15          | 170                | 108                | 62           | 36%                    | 57%                    | 64%                          |
| 16          | 510                | 381                | 129          | 25%                    | 34%                    | 75%                          |
| 17          | 200                | 67                 | 133          | 67%                    | 199%                   | 34%                          |
| 18          | 170                | 109                | 61           | 36%                    | 56%                    | 64%                          |
| 19          | 360                | 199                | 161          | 45%                    | 81%                    | 55%                          |
| 20          | 260                | 164                | 96           | 37%                    | 59%                    | 63%                          |
| 21          | 440                | 372                | 68           | 15%                    | 18%                    | 85%                          |

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN<br>to<br>Total PTH | PIN<br>to<br>Whole PTH | Whole PTH<br>to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|------------------------|------------------------|------------------------------|
| 22          | 120                | 51.7               | 68.3         | 57%                    | 132%                   | 43%                          |
| 23          | 600                | 527                | 73           | 12%                    | 14%                    | 83%                          |
| 24          | 220                | 130                | 90           | 41%                    | 69%                    | 59%                          |
| 25          | 190                | 136                | 54           | 28%                    | 40%                    | 72%                          |
| 26          | 220                | 118                | 102          | 46%                    | 86%                    | 54%                          |
| 27          | 630                | 334                | 296          | 47%                    | 89%                    | 53%                          |
| 28          | 150                | 90                 | 60           | 40%                    | 67%                    | 60%                          |
| 29          | 170                | 106                | 64           | 38%                    | 60%                    | 62%                          |
| 30          | 810                | 489                | 321          | 40%                    | 66%                    | 60%                          |
| 31          | 570                | 319                | 251          | 44%                    | 79%                    | 56%                          |
| 32          | 570                | 467                | 103          | 18%                    | 22%                    | 82%                          |
| 33          | 400                | 300                | 100          | 25%                    | 33%                    | 75%                          |
| 34          | 560                | 378                | 182          | 33%                    | 48%                    | 68%                          |
| 35          | 310                | 121                | 189          | 61%                    | 156%                   | 39%                          |
| 36          | 240                | 98                 | 142          | 59%                    | 145%                   | 41%                          |
| 37          | 280                | 133                | 157          | 54%                    | 118%                   | 48%                          |
| 38          | 230                | 124                | 106          | 46%                    | 85%                    | 54%                          |
| 39          | 350                | 319                | 31           | 9%                     | 10%                    | 91%                          |
| 40          | 200                | 133                | 67           | 34%                    | 50%                    | 67%                          |
| 41          | 920                | 564                | 356          | 39%                    | 63%                    | 61%                          |
| 42          | 210                | 89                 | 121          | 58%                    | 136%                   | 42%                          |
| 43          | 1990               | 904                | 1086         | 55%                    | 120%                   | 45%                          |
| 44          | 300                | 212                | 88           | 29%                    | 42%                    | 71%                          |
| 45          | 260                | 132                | 128          | 49%                    | 97%                    | 51%                          |
| 46          | 140                | 72                 | 68           | 49%                    | 94%                    | 51%                          |
| 47          | 250                | 129                | 121          | 48%                    | 94%                    | 52%                          |
| 48          | 130                | 72                 | 58           | 45%                    | 81%                    | 56%                          |
| 49          | 1840               | 1000               | 840          | 46%                    | 84%                    | 54%                          |
| 50          | 280                | 167                | 113          | 40%                    | 68%                    | 60%                          |
| 51          | 490                | 268                | 222          | 45%                    | 83%                    | 55%                          |
| 52          | 150                | 77.1               | 72.9         | 49%                    | 95%                    | 51%                          |
| 53          | 140                | 58.1               | 81.9         | 59%                    | 141%                   | 42%                          |
| 54          | 210                | 92.7               | 117.3        | 56%                    | 127%                   | 44%                          |
| 55          | 160                | 79                 | 81           | 51%                    | 103%                   | 49%                          |
| 56          | 480                | 296                | 184          | 38%                    | 62%                    | 62%                          |
| 57          | 480                | 281                | 199          | 41%                    | 71%                    | 59%                          |
| 58          | 270                | 120                | 150          | 56%                    | 125%                   | 44%                          |

5

10

15

20

25

30

35

40

45

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN<br>to<br>Total PTH | PIN<br>to<br>Whole PTH | Whole PTH<br>to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|------------------------|------------------------|------------------------------|
| 59          | 97                 | 45                 | 52           | 54%                    | 116%                   | 46%                          |
| 60          | 330                | 154                | 176          | 53%                    | 114%                   | 47%                          |
| 61          | 110                | 56                 | 54           | 49%                    | 96%                    | 51%                          |
| 62          | 660                | 456                | 204          | 31%                    | 45%                    | 69%                          |
| 63          | 300                | 137                | 163          | 54%                    | 119%                   | 46%                          |
| 64          | 240                | 145                | 95           | 40%                    | 66%                    | 60%                          |
| 65          | 100                | 66.5               | 33.5         | 34%                    | 50%                    | 67%                          |
| 66          | 410                | 416.3              | -6.3         | -2%                    | -2%                    | 102%                         |
| 67          | 410                | 235.7              | 174.3        | 43%                    | 74%                    | 57%                          |
| 68          | 45                 | 14.4               | 30.6         | 68%                    | 213%                   | 32%                          |
| 69          | 200                | 102.3              | 97.7         | 49%                    | 96%                    | 51%                          |
| 70          | 300                | 134                | 166          | 55%                    | 124%                   | 45%                          |
| 71          | 320                | 202                | 118          | 37%                    | 58%                    | 63%                          |
| 72          | 440                | 254                | 186          | 42%                    | 73%                    | 58%                          |
| 73          | 190                | 99.6               | 90.4         | 48%                    | 91%                    | 52%                          |
| 74          | 160                | 74.6               | 85.4         | 53%                    | 114%                   | 47%                          |
| 75          | 600                | 429.8              | 170.2        | 28%                    | 40%                    | 72%                          |
| 76          | 1140               | 632                | 508          | 45%                    | 80%                    | 55%                          |
| 77          | 440                | 211                | 229          | 52%                    | 109%                   | 48%                          |
| 78          | 450                | 276                | 174          | 39%                    | 63%                    | 61%                          |
| 79          | 510                | 344                | 166          | 33%                    | 48%                    | 67%                          |
| 80          | 190                | 62.8               | 127.2        | 67%                    | 203%                   | 33%                          |
| 81          | 170                | 86                 | 84           | 49%                    | 98%                    | 51%                          |
| 82          | 180                | 103.4              | 76.6         | 43%                    | 74%                    | 57%                          |
| 83          | 78                 | 22.7               | 55.3         | 71%                    | 244%                   | 29%                          |
| 84          | 230                | 117                | 113          | 49%                    | 97%                    | 51%                          |
| 85          | 160                | 96                 | 64           | 40%                    | 67%                    | 60%                          |
| 86          | 220                | 89                 | 131          | 60%                    | 147%                   | 40%                          |
| 87          | 470                | 321.5              | 148.5        | 32%                    | 46%                    | 68%                          |
| 88          | 310                | 137                | 173          | 56%                    | 126%                   | 44%                          |
| 89          | 2050               | 1127               | 923          | 45%                    | 82%                    | 55%                          |
| 90          | 930                | 414                | 516          | 55%                    | 125%                   | 45%                          |
| 91          | 180                | 65                 | 115          | 64%                    | 177%                   | 36%                          |
| 92          | 560                | 238                | 322          | 58%                    | 135%                   | 43%                          |
| 93          | 640                | 597                | 43           | 7%                     | 7%                     | 93%                          |
| 94          | 590                | 382                | 208          | 35%                    | 54%                    | 65%                          |
| 95          | 270                | 103                | 167          | 62%                    | 162%                   | 38%                          |

50

55

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN<br>to<br>Total PTH | PIN<br>to<br>Whole PTH | Whole PTH<br>to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|------------------------|------------------------|------------------------------|
| 96          | 560                | 349                | 211          | 38%                    | 60%                    | 62%                          |
| 97          | 180                | 78                 | 102          | 57%                    | 131%                   | 43%                          |
| 98          | 790                | 429                | 361          | 46%                    | 84%                    | 54%                          |
| 99          | 670                | 372                | 298          | 44%                    | 80%                    | 56%                          |
| 100         | 140                | 20.4               | 119.6        | 85%                    | 586%                   | 15%                          |
| 101         | 190                | 117                | 73           | 38%                    | 62%                    | 62%                          |
| 102         | 190                | 108                | 82           | 43%                    | 76%                    | 57%                          |
| 103         | 430                | 217                | 213          | 50%                    | 98%                    | 50%                          |
| 104         | 560                | 439                | 121          | 22%                    | 28%                    | 78%                          |
| 105         | 500                | 357.7              | 142.3        | 28%                    | 40%                    | 72%                          |
| 106         | 1560               | 777                | 783          | 50%                    | 101%                   | 50%                          |
| 107         | 62                 | 24.3               | 37.7         | 61%                    | 155%                   | 39%                          |
| 108         | 430                | 226                | 204          | 47%                    | 90%                    | 53%                          |
| 109         | 160                | 67.2               | 92.8         | 58%                    | 138%                   | 42%                          |
| 110         | 530                | 346                | 184          | 35%                    | 53%                    | 65%                          |
| 111         | 260                | 142                | 118          | 45%                    | 83%                    | 55%                          |
| 112         | 580                | 163                | 417          | 72%                    | 256%                   | 28%                          |
| 113         | 440                | 579                | -139         | -32%                   | -24%                   | 132%                         |
| 114         | 500                | 232.3              | 267.7        | 54%                    | 115%                   | 46%                          |
| 115         | 160                | 60                 | 100          | 63%                    | 167%                   | 38%                          |
| 116         | 340                | 202                | 138          | 41%                    | 68%                    | 59%                          |
| 117         | 260                | 138                | 122          | 47%                    | 88%                    | 53%                          |
| 118         | 260                | 119                | 141          | 54%                    | 118%                   | 46%                          |
| 119         | 160                | 84                 | 76           | 48%                    | 90%                    | 53%                          |
| 120         | 130                | 46                 | 84           | 65%                    | 183%                   | 35%                          |
| 121         | 190                | 104                | 86           | 45%                    | 83%                    | 55%                          |
| 122         | 420                | 334                | 86           | 20%                    | 26%                    | 80%                          |
| 123         | 630                | 440                | 190          | 30%                    | 43%                    | 70%                          |
| 124         | 75                 | 26.4               | 48.6         | 65%                    | 184%                   | 35%                          |
| 125         | 260                | 143                | 117          | 45%                    | 82%                    | 55%                          |
| 126         | 640                | 409                | 231          | 36%                    | 56%                    | 64%                          |
| 127         | 130                | 66.7               | 63.3         | 49%                    | 95%                    | 51%                          |
| 128         | 700                | 381                | 319          | 46%                    | 84%                    | 54%                          |
| 129         | 560                | 376                | 184          | 33%                    | 49%                    | 67%                          |
| 130         | 240                | 107                | 133          | 55%                    | 124%                   | 45%                          |
| 131         | 110                | 63                 | 47           | 43%                    | 75%                    | 57%                          |
| 132         | 420                | 297                | 123          | 29%                    | 41%                    | 71%                          |

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN to<br>Total PTH | PIN to<br>Whole PTH | Whole PTH to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|---------------------|---------------------|---------------------------|
| 133         | 580                | 229                | 351          | 61%                 | 153%                | 39%                       |
| 134         | 310                | 201.2              | 108.8        | 35%                 | 54%                 | 65%                       |
| 135         | 160                | 97.9               | 62.1         | 39%                 | 63%                 | 61%                       |
| 136         | 290                | 138.7              | 151.3        | 52%                 | 109%                | 48%                       |
| 137         | 200                | 96.2               | 103.8        | 52%                 | 108%                | 48%                       |
| 138         | 770                | 662.7              | 107.3        | 14%                 | 16%                 | 86%                       |
| 139         | 290                | 130.7              | 159.3        | 55%                 | 122%                | 45%                       |
| 140         | 260                | 219                | 41           | 16%                 | 19%                 | 84%                       |
| 141         | 350                | 211                | 139          | 40%                 | 66%                 | 60%                       |
| 142         | 730                | 463.5              | 266.5        | 37%                 | 57%                 | 63%                       |
| 143         | 490                | 231                | 259          | 53%                 | 112%                | 47%                       |
| 144         | 160                | 87                 | 73           | 46%                 | 84%                 | 54%                       |
| 145         | 380                | 222                | 158          | 42%                 | 71%                 | 58%                       |
| 146         | 210                | 93.5               | 116.5        | 55%                 | 125%                | 45%                       |
| 147         | 630                | 383.4              | 246.6        | 39%                 | 64%                 | 61%                       |
| 148         | 150                | 83.2               | 66.8         | 45%                 | 80%                 | 55%                       |
| 149         | 320                | 152.5              | 167.5        | 52%                 | 110%                | 48%                       |
| 150         | 900                | 467.6              | 432.4        | 48%                 | 92%                 | 52%                       |
| 151         | 1180               | 818.6              | 361.4        | 31%                 | 44%                 | 69%                       |
| 152         | 120                | 38.4               | 81.6         | 68%                 | 213%                | 32%                       |
| 153         | 5230               | 1388               | 3842         | 73%                 | 277%                | 27%                       |
| 154         | 34                 | 10.5               | 23.5         | 69%                 | 224%                | 31%                       |
| 155         | 1020               | 590.6              | 429.4        | 42%                 | 73%                 | 58%                       |
| 156         | 180                | 76.6               | 103.4        | 57%                 | 135%                | 43%                       |
| 157         | 120                | 51.1               | 68.9         | 57%                 | 135%                | 43%                       |
| Median      | 300                | 154                | 127          | 46%                 | 84%                 | 54%                       |

TABLE 2 shows the results, individually and comparatively, of the wPTH, PIN, and total PTH assays from the normals.

TABLE 2

| Patient No. | Total PTH<br>pg/ml | Whole PTH<br>pg/ml | PIN<br>pg/ml | PIN to<br>Total PTH | PIN to<br>Whole PTH | Whole PTH to<br>Total PTH |
|-------------|--------------------|--------------------|--------------|---------------------|---------------------|---------------------------|
| 1           | 17.13              | 3.32               | 13.81        | 81%                 | 416%                | 19%                       |
| 2           | 32.92              | 10.49              | 22.43        | 68%                 | 214%                | 32%                       |
| 3           | 31.32              | 10.31              | 21.01        | 67%                 | 204%                | 33%                       |
| 4           | 41.84              | 12.72              | 29.12        | 70%                 | 229%                | 30%                       |
| 5           | 33.03              | 10.09              | 22.94        | 69%                 | 227%                | 31%                       |
| 6           | 44.32              | 14.23              | 30.09        | 68%                 | 211%                | 32%                       |
| 7           | 31.47              | 6.8                | 24.67        | 78%                 | 363%                | 22%                       |
| 8           | 20.82              | 10.03              | 10.79        | 52%                 | 108%                | 48%                       |
| 9           | 34.64              | 15.95              | 18.69        | 54%                 | 117%                | 46%                       |
| 10          | 23.69              | 5.25               | 18.44        | 78%                 | 351%                | 22%                       |
| 11          | 53.98              | 17.82              | 36.16        | 67%                 | 203%                | 33%                       |
| 12          | 52.71              | 18.83              | 33.88        | 64%                 | 180%                | 36%                       |
| 13          | 26.92              | 5.63               | 21.29        | 79%                 | 378%                | 21%                       |
| 14          | 39.93              | 11.86              | 28.07        | 70%                 | 237%                | 30%                       |
| 15          | 48.84              | 20.47              | 28.37        | 58%                 | 139%                | 42%                       |
| 16          | 29.56              | 13.68              | 15.88        | 54%                 | 116%                | 46%                       |
| 17          | 36.19              | 14.69              | 21.5         | 59%                 | 146%                | 41%                       |
| 18          | 20.96              | 6.99               | 13.97        | 67%                 | 200%                | 33%                       |
| 19          | 59.29              | 27.89              | 31.4         | 53%                 | 113%                | 47%                       |
| 20          | 45.57              | 18.23              | 27.34        | 60%                 | 150%                | 40%                       |
| 21          | 35.64              | 18.72              | 16.92        | 47%                 | 90%                 | 53%                       |
| 22          | 38.53              | 19.56              | 18.97        | 49%                 | 97%                 | 51%                       |
| 23          | 21.71              | 9.34               | 12.37        | 57%                 | 132%                | 43%                       |
| 24          | 32.42              | 13.51              | 18.91        | 58%                 | 140%                | 42%                       |
| 25          | 28.5               | 10.41              | 18.09        | 63%                 | 174%                | 37%                       |
| 26          | 18.17              | 7.8                | 10.37        | 57%                 | 133%                | 43%                       |
| 27          | 39.96              | 17.29              | 22.67        | 57%                 | 131%                | 43%                       |
| 28          | 34.08              | 15.24              | 18.84        | 55%                 | 124%                | 45%                       |
| 29          | 42.95              | 19.59              | 23.36        | 54%                 | 119%                | 46%                       |
| 30          | 38.4               | 12.16              | 26.24        | 68%                 | 216%                | 32%                       |
| 31          | 47.57              | 18.45              | 29.12        | 61%                 | 158%                | 39%                       |
| Median      | 34.64              | 13.51              | 21.5         | 61%                 | 158%                | 39%                       |



Clearly, the statistically significant differences in the medians of these two groups demonstrates that one can differentiate between the two by using these assays alone or by comparing their respective values.

TABLE 3

| Sample Type                    | Total PTH (pg/mL) | Whole PTH (pg/mL) | PIN (pg/mL) | PIN to Total PTH | PIN to Whole PTH | Whole PTH to Total PTH |
|--------------------------------|-------------------|-------------------|-------------|------------------|------------------|------------------------|
| Chronic Uremia (n=157) Medians | 300               | 154               | 127         | 46%              | 84%              | 55%                    |
| Normal (n=31) Medians          | 34.64             | 13.51             | 21.5        | 61%              | 158%             | 37%                    |
| P-Value                        | < 0.0001          | < 0.0001          | < 0.0001    | < 0.0001         | < 0.0001         | < 0.0001               |

The ordinarily skilled artisan can appreciate that the present invention can incorporate any number of the preferred features described above.

All publications or unpublished patent applications mentioned herein are hereby incorporated by reference thereto.

Other embodiments of the present invention are not presented here which are obvious to those of ordinary skill in the art, now or during the term of any patent issuing from this patent specification, and thus, are within the spirit and scope of the present invention.

## Claims

5

10

15

20

25

30

35

40

45

50

55

**WE CLAIM:**

1. A method for differentiating between a person having substantially normal parathyroid function and having hyperparathyroidism comprising determining and comparing at least two of the parameters selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and the total parathyroid hormone level in the person.
2. The method of Claim 1 wherein the comparison is in the form of a ratio or proportion.
3. The method of Claim 1 wherein the person is a patient with chronic uremia.
4. The method of Claim 1 wherein one measures the whole parathyroid hormone level and the total parathyroid hormone level in the person, determines the parathyroid hormone inhibitory peptide fragment level from these two measurements, and compares the whole parathyroid hormone level to the parathyroid hormone inhibitory peptide fragment level.
5. The method of Claim 4 wherein the comparison is in the form of a ratio or proportion.
6. The method of Claim 4 wherein the person is a patient with chronic uremia.
7. A method for differentiating between a person having substantially normal parathyroid function and having hyperparathyroidism comprising determining one parameter selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and a calculated total parathyroid hormone level.

5 8. The method of Claim 7 wherein one determines the parathyroid hormone inhibitory peptide fragment level by measuring the whole parathyroid hormone level and the total parathyroid hormone level.

10 9. The method of Claim 7 wherein one determines the total parathyroid hormone level by measuring the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.

15 10. The method of Claim 7 wherein the person is a patient with chronic uremia.

10 11. The method of Claim 1 wherein one measures and compares the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.

20 12. The method of Claim 11 wherein the comparison is in the form of a ratio or proportion.

25 13. The method of Claim 11 wherein the person is a patient with chronic uremia.

30 14. The method of Claim 1 wherein one measures and compares the whole parathyroid hormone level and the total parathyroid hormone level in the person.

35 15. The method of Claim 14 wherein the comparison is in the form of a ratio or proportion.

40 16. The method of Claim 14 wherein the person is a patient with chronic uremia.

45 17. The method of Claim 1 wherein one measures and compares the parathyroid hormone inhibitory peptide fragment level and the total parathyroid hormone level in the person.

5 18. The method of Claim 17 wherein the comparison is in the form of a ratio or proportion.

10 19. The method of Claim 17 wherein the person is a patient with chronic uremia.

5 20. A method for monitoring parathyroid related bone diseases and treatments therefor comprising determining and comparing at least two of the parameters selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and the total parathyroid hormone level in the person.

15 21. The method of Claim 20 wherein the comparison is in the form of a ratio or proportion.

20 22. The method of Claim 20 wherein one measures the whole parathyroid hormone level and the total parathyroid hormone level in the person, determines the parathyroid hormone inhibitory peptide fragment level from these two measurements, and compares the whole parathyroid hormone level to the parathyroid hormone inhibitory peptide fragment level.

25 23. The method of Claim 22 wherein the comparison is in the form of a ratio or proportion.

30 24. The method of Claim 20 wherein one measures and compares the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.

35 25. The method of Claim 20 wherein one measures and compares the whole parathyroid hormone level and the total parathyroid hormone level in the person.

- 5 26. The method of Claim 24 wherein the comparison is in the form of a ratio or proportion.
- 10 27. The method of Claim 25 wherein the comparison is in the form of a ratio or proportion.
- 15 28. The method of Claim 20 wherein one measures and compares the parathyroid hormone inhibitory peptide fragment level and the total parathyroid hormone level in the person.
- 20 29. The method of Claim 28 wherein one measures the whole parathyroid hormone level in order to calculate the parathyroid hormone inhibitory peptide fragment level from the whole parathyroid hormone level and the total parathyroid hormone level.
- 25 30. The method of Claim 28 wherein the comparison is in the form of a ratio or proportion.
- 30 31. A method for monitoring parathyroid related bone diseases and treatments therefor comprising determining one parameter selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and the calculated total parathyroid hormone level.
- 35 32. The method of Claim 31 wherein one determines the parathyroid hormone inhibitory peptide fragment level by measuring the whole parathyroid hormone level and the total parathyroid hormone level.
- 40 33. The method of Claim 31 wherein one determines the total parathyroid hormone level by measuring the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.
- 45

5  
10  
15  
20  
25  
30  
35  
40  
45  
50  
55

34. A method for monitoring the effects of therapeutic treatment for hyperparathyroidism comprising determining and comparing at least two of the parameters selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and the total parathyroid hormone level in the person.

35. The method of Claim 34 wherein the comparison is in the form of a ratio or proportion.

36. The method of Claim 34 wherein one measures the whole parathyroid hormone level and the total parathyroid hormone level in the person, determines the parathyroid hormone inhibitory peptide fragment level from these two measurements, and compares the whole parathyroid hormone level to the parathyroid hormone inhibitory peptide fragment level.

37. The method of Claim 36 wherein the comparison is in the form of a ratio or proportion.

38. The method of Claim 34 wherein one measures and compares the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.

39. The method of Claim 34 wherein one measures and compares the whole parathyroid hormone level and the total parathyroid hormone level in the person.

40. The method of Claim 39 wherein the comparison is in the form of a ratio or proportion.

41. The method of Claim 34 wherein one determines and compares the parathyroid

- 5 hormone inhibitory peptide fragment level and the total parathyroid hormone level in the person.
- 10 42. The method of Claim 41 wherein one determines the parathyroid hormone inhibitory fragment level by measuring the whole parathyroid hormone and the total parathyroid hormone level.
- 15 43. The method of Claim 41 wherein the comparison is in the form of a ratio or proportion.
- 10 44. A method for monitoring the effects of therapeutic treatment for hyperparathyroidism comprising determining one parameter selected from the group consisting of the whole parathyroid hormone level, the parathyroid hormone inhibitory peptide fragment level, and the calculated total parathyroid hormone level.
- 20 45. The method of Claim 44 wherein one determines the parathyroid hormone inhibitory peptide fragment level by measuring the whole parathyroid hormone level and the total parathyroid hormone level.
- 25 46. The method of Claim 44 wherein one determines the total parathyroid hormone level by measuring the whole parathyroid hormone level and the parathyroid hormone inhibitory peptide fragment level.



### Whole Human PTH (1-84)

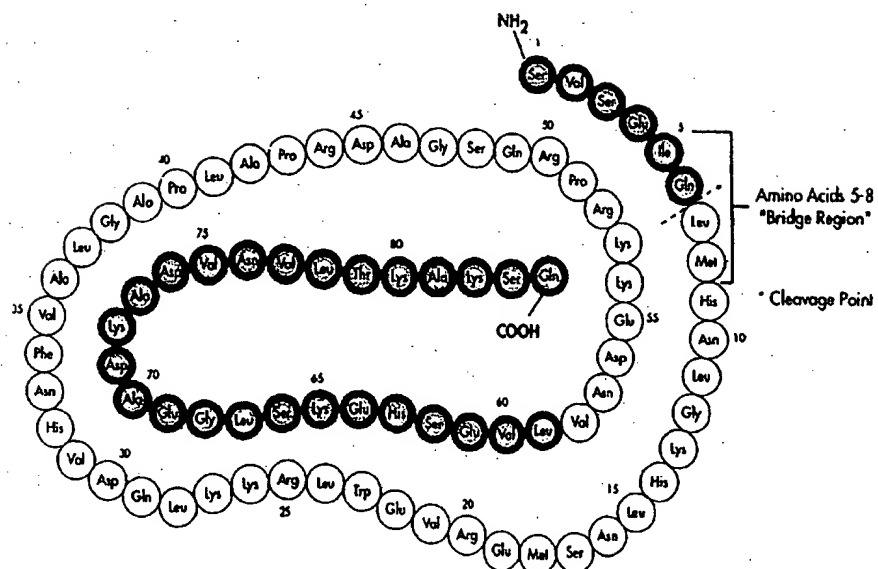


FIG. 2



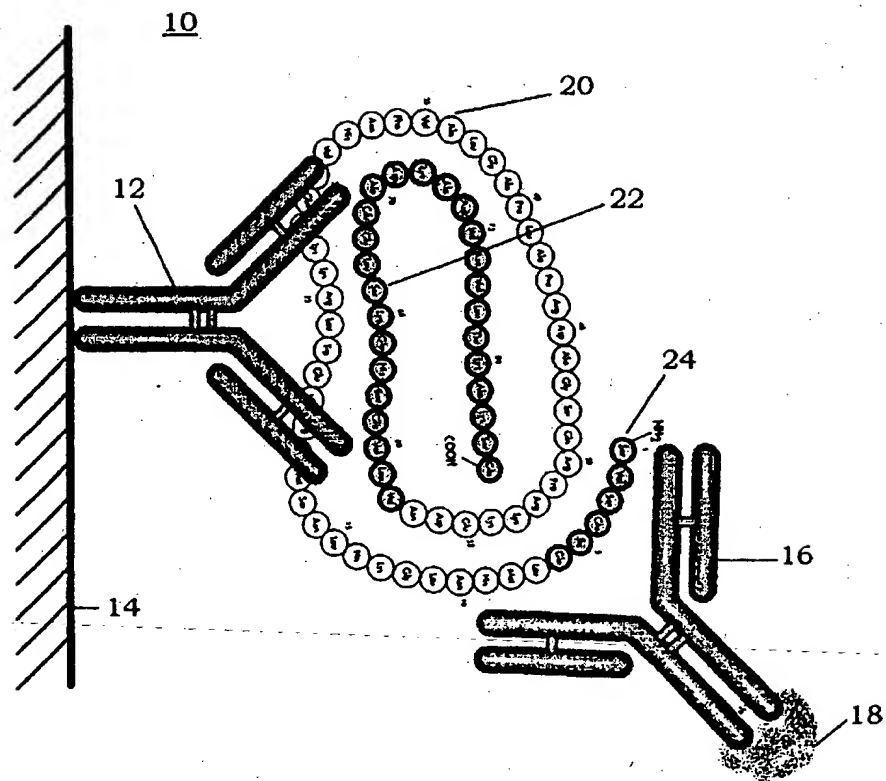


FIG. 3

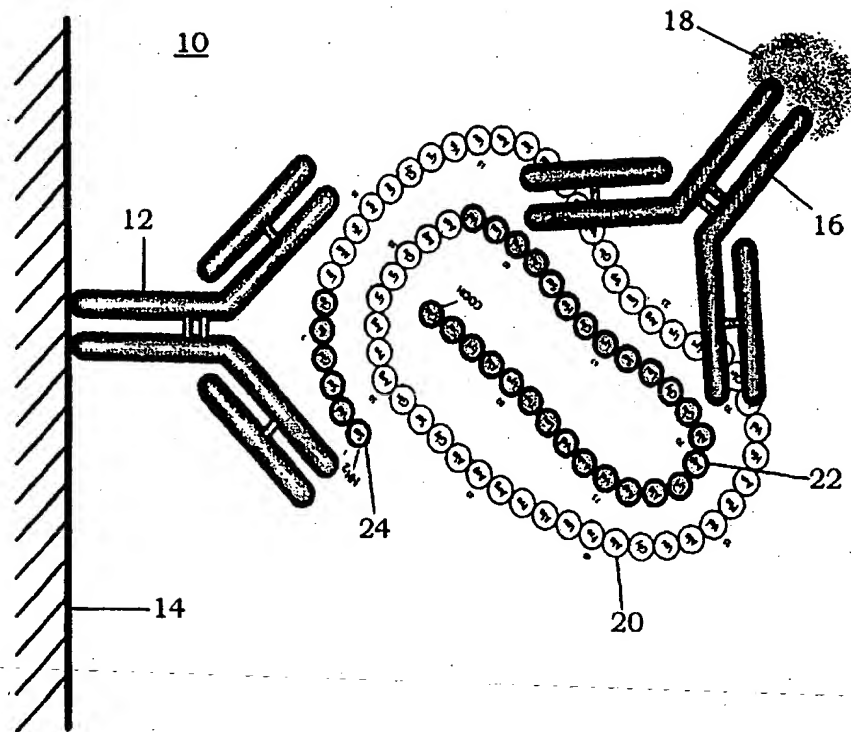

# FIG. 4

## Standard Curve for Whole PTH Assay

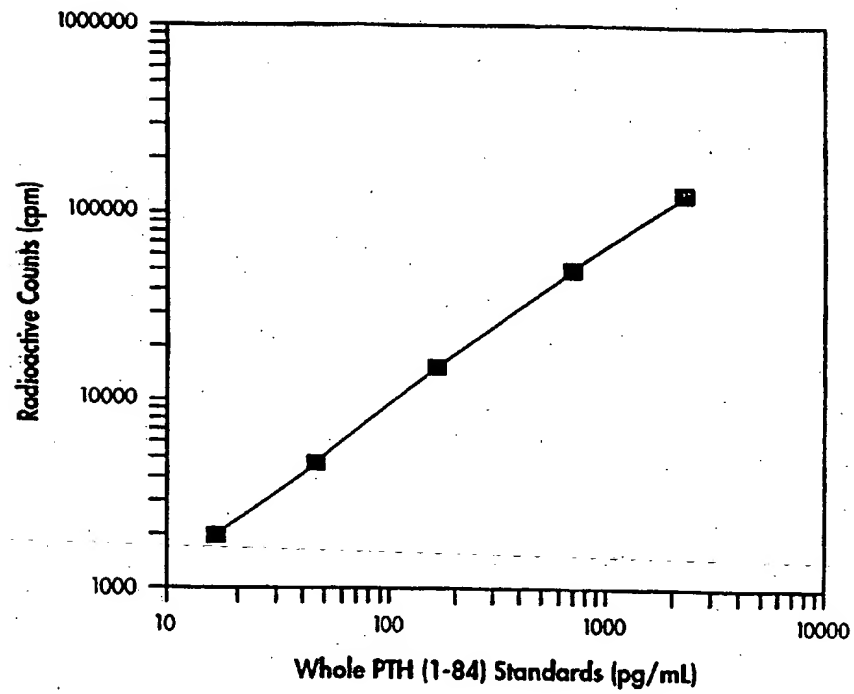

## FIG. 5

**Normal Value Comparison**  
Whole PTH Assay (with PTH 1-8 Antibody as Tracer)  
versus  
Nichols' Intact PTH Assay (with PTH 7-84 Interference)

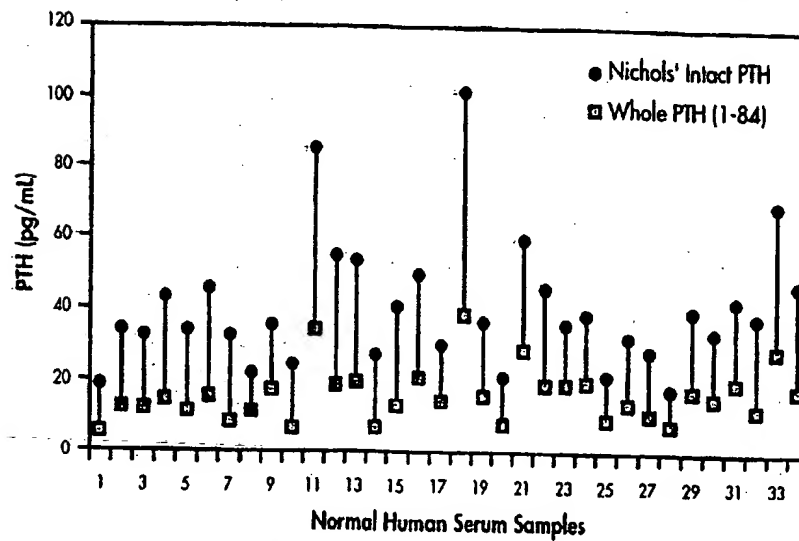


FIG. 6

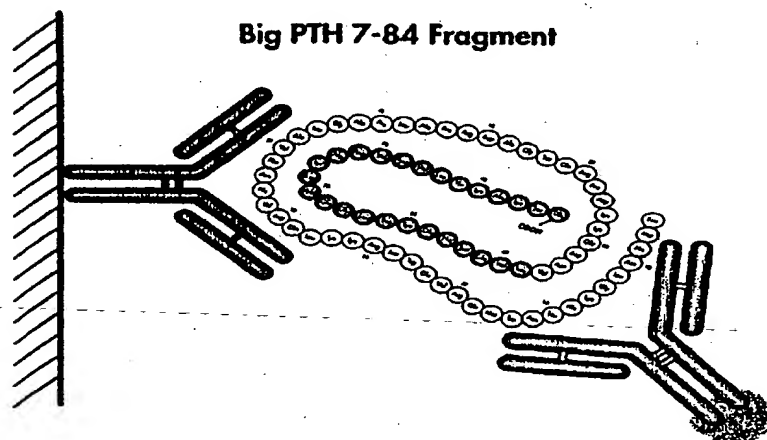
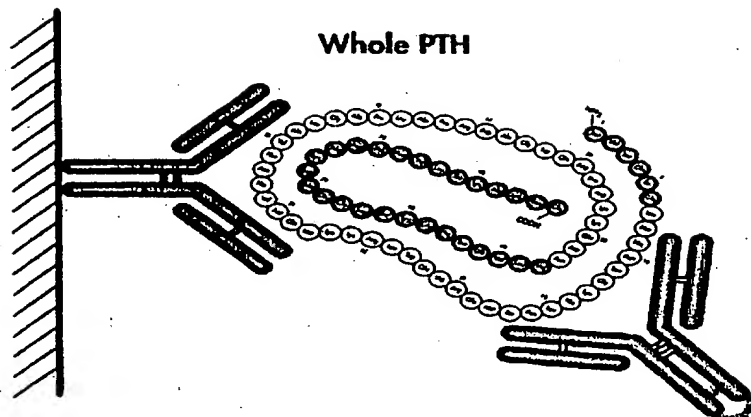


FIG. 7

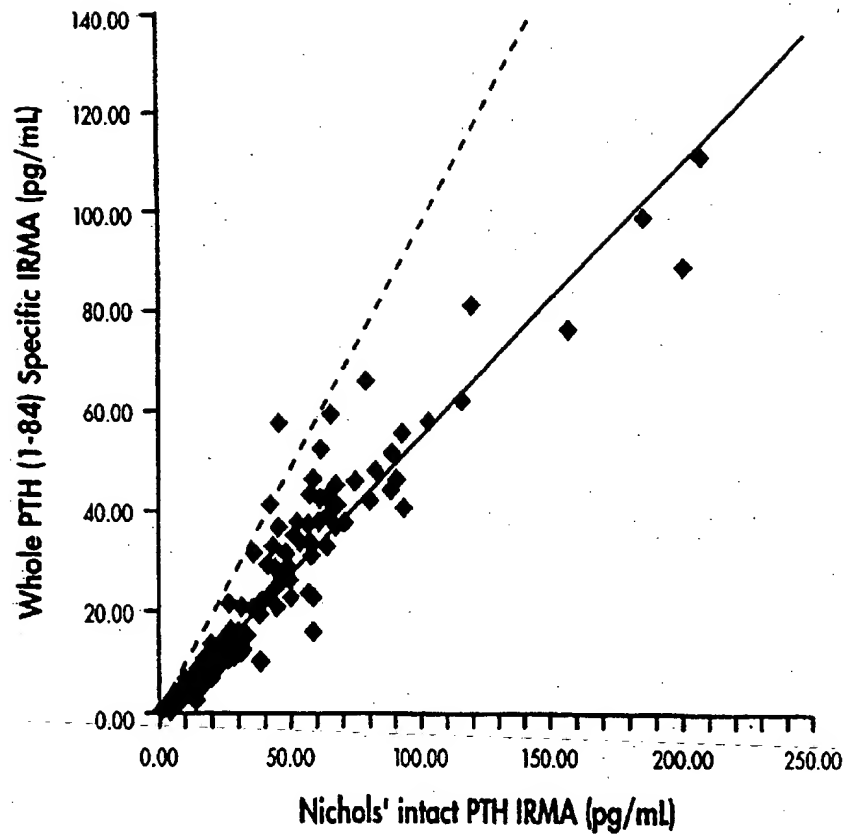


FIG. 8

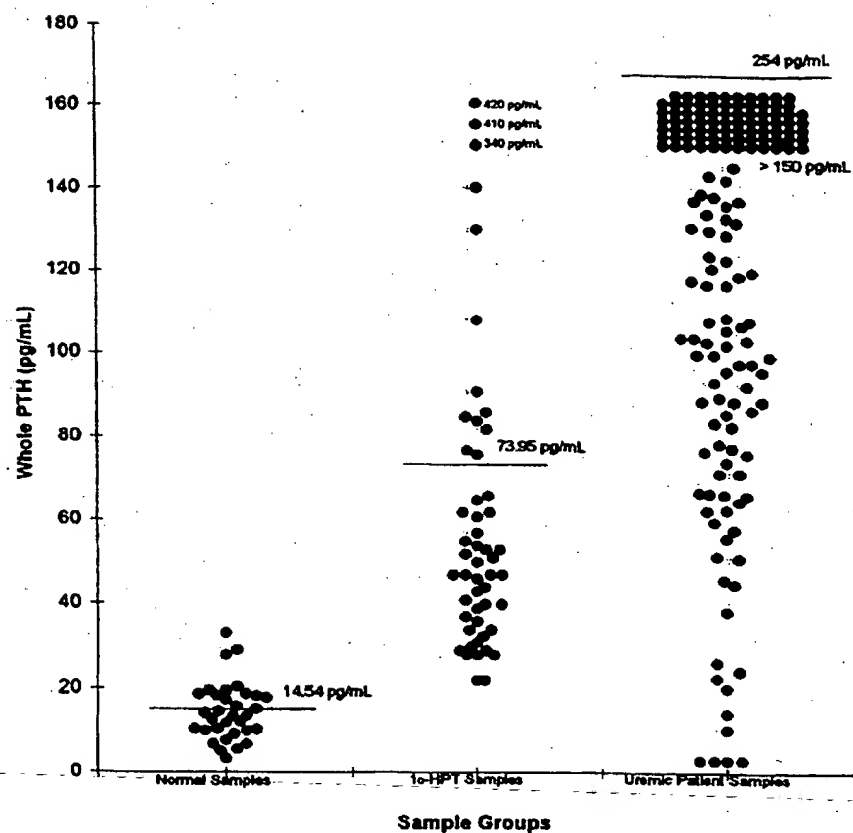




FIG. 9

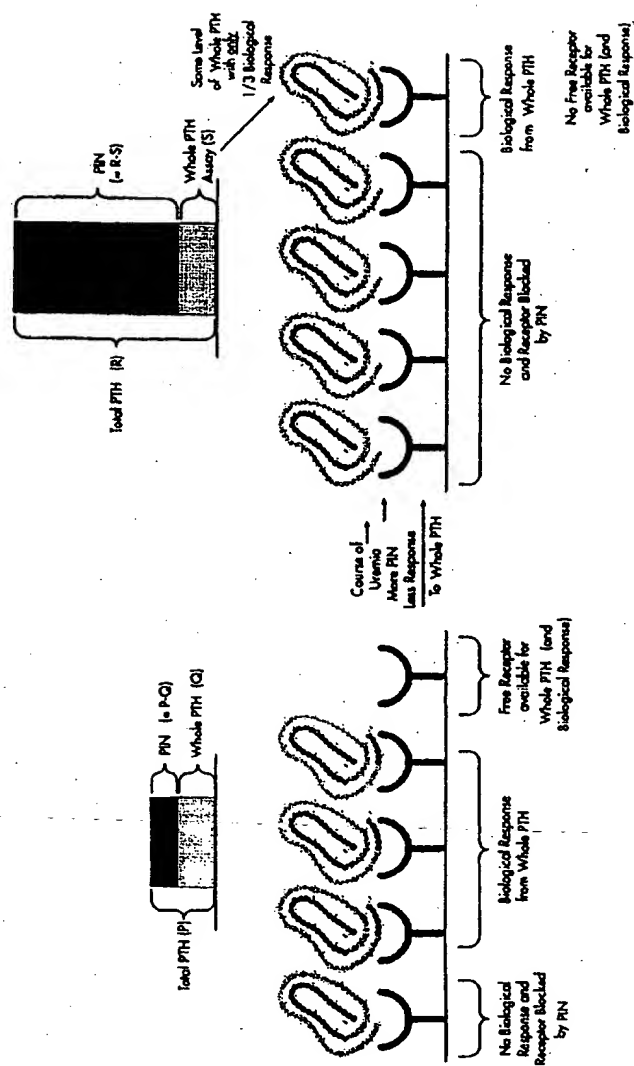


FIG. 10

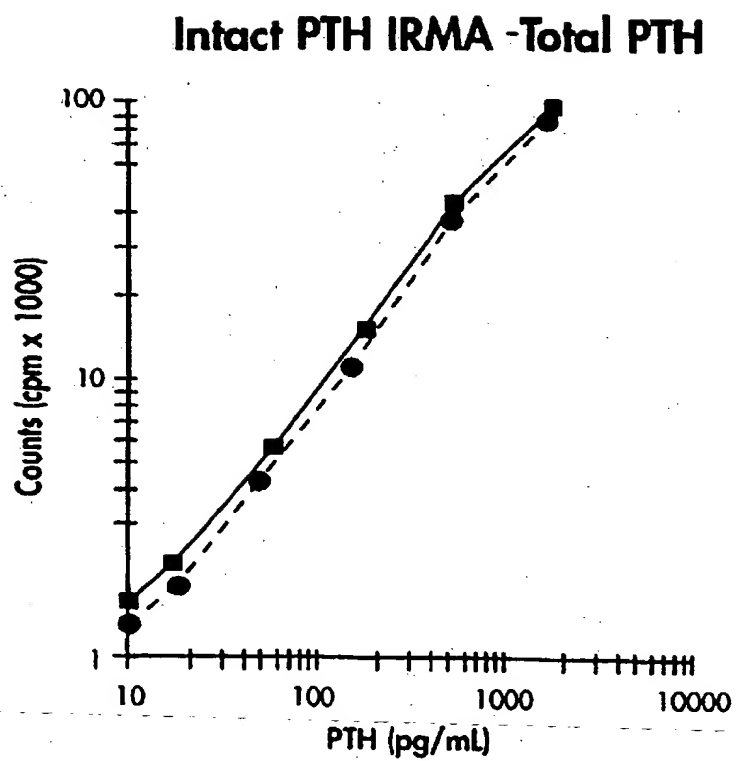


FIG. 11


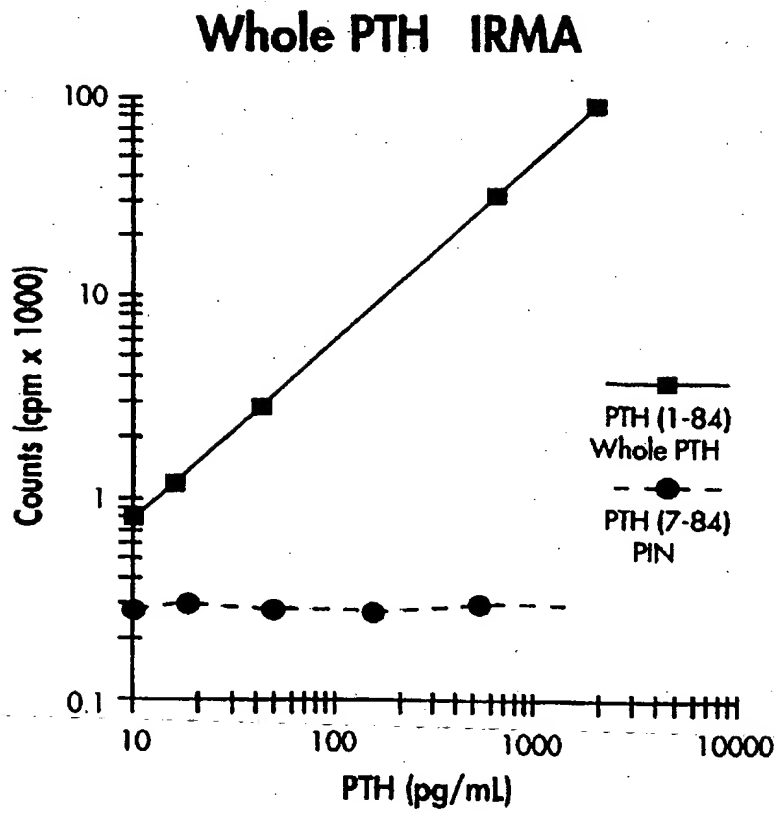


FIG. 12

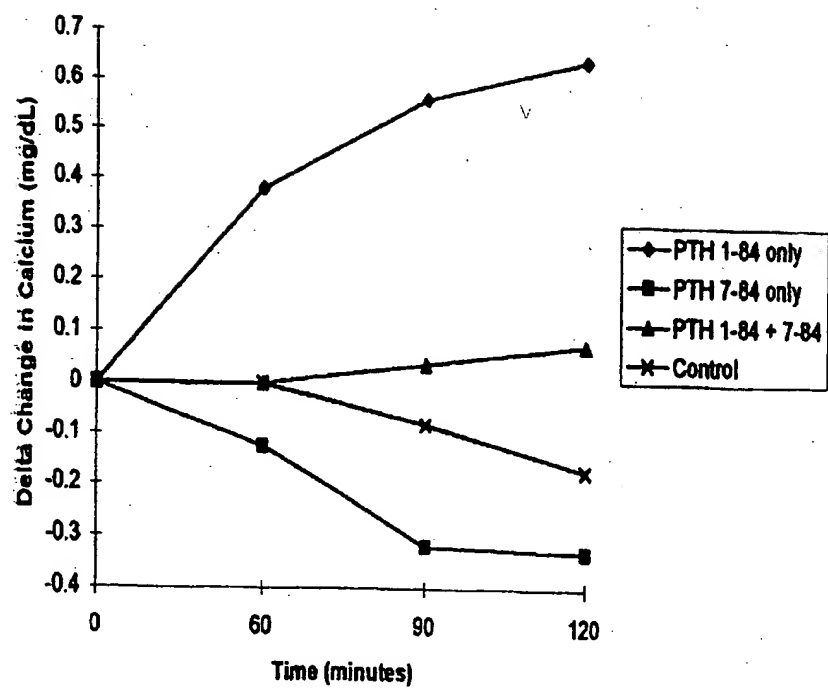

## SEQUENCE LISTING

<110> Cantor, Thomas L.  
Gao, Ping

<120> Methods for Differentiating Parathyroid and Bone Status Related Diseases

<160> 5

<170> Microsoft Word 7.0

<210> 1

<211> 84 [integer length]

<212> PRT

<400> 1

Ser Val Ser Glu Ile Gln Leu Met His Asn Leu Gly Lys His Leu

1 5 10 15

Asn Ser Met Glu Arg Val Glu Trp Leu Arg Lys Lys Leu Gln Asp

20 25 30

Val His Asn Phe Val Ala Leu Gly Ala Pro Leu Ala Pro Arg Asp

35 40 45

Ala Gly Ser Gln Arg Pro Arg Lys Lys Glu Asp Asn Val Leu Val

50 55 60

Glu Ser His Glu Lys Ser Leu Gly Glu Ala Asn Lys Ala Asp Val  
65 70 75

Asn Val Leu Thr Lys Ala Lys Ser Gln  
80

<210> 2

<211> 82 [integer length]

<212> PRT

<400> 2

Ser Glu Ile Gln Leu Met His Asn Leu Gly Lys His Leu Asn Ser  
1 5 10 15

Met Glu Arg Val Glu Trp Leu Arg Lys Lys Leu Gln Asp Val His  
20 25 30

Asn Phe Val Ala Leu Gly Ala Pro Leu Ala Pro Arg Asp Ala Gly  
35 40 45

Ser Gln Arg Pro Arg Lys Lys Glu Asp Asn Val Leu Val Glu Ser  
50 55 60

His Glu Lys Ser Leu Gly Glu Ala Asn Lys Ala Asp Val Asn Val  
65 70 75

Leu Thr Lys Ala Lys Ser Gln  
80

<210> 3

<211> 51 [integer length]

<212> PRT

<400> 3.

Phe Val Ala Leu Gly Ala Pro Leu Ala Pro Arg Asp Ala Gly Ser  
1 5 10 15

Gln Arg Pro Arg Lys Lys Glu Asp Asn Val Leu Val Glu Ser His  
20 25 30

Glu Lys Ser Leu Gly Glu Ala Asn Lys Ala Asp Val Asn Val Leu  
35 40 45

Thyr Lys Ala Lys Ser Gln  
50

<210> 4

<211> 8 [integer length]

<212> PRT

<400> 4

Ser Val Ser Glu Ile Gln Leu Met

1 5

<210> 5

<211> 32 [integer length]

<212> PRT

<400> 5

Leu Met His Asn Leu Gly Lys His Leu Asn Ser Met Glu Arg Val

1 5 10 15

Glu Trp Leu Arg Lys Lys Leu Gln Asp Val His Asn Phe Val Ala

20 25 30

Leu Gly



## INTERNATIONAL SEARCH REPORT

International application No.  
PCT/US00/00855

## A. CLASSIFICATION OF SUBJECT MATTER

IPC(7) : G01N 33/74

US CL : 436/87

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 436/87, 518, 536, 548, 811; 435/7.94; 530/388.24, 389.2

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

DIALOG

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

| Category* | Citation of document, with indication, where appropriate, of the relevant passages   | Relevant to claim No. |
|-----------|--|-----------------------|
| X         | BROSSARD et al. Accumulation of a Non-(1-84) Molecular Form of Parathyroid Hormone (PTH) Detected by Intact PTH Assay in Renal Failure: Importance in the Interpretation of PTH Values. Journal of Clinical Endocrinology and Metabolism. 1996, Vol. 81, No. 11, pages 3923-3929, see entire document. | 1-46                  |
| Y         | LEPAGE et al. A Non-(1-84) Circulating Parathyroid Hormone (PTH) Fragment Interferes Significantly with Intact PTH Commercial Assay Measurements in Uremic Samples. Clinical Chemistry. April 1998, Vol. 44, No. 4, pages 805-809; see entire document.  | 1-46                  |

☒ Further documents are listed in the continuation of Box C. ☐ See patent family annex.

|   |  |
|---|--|
| * Special categories of cited documents:  | *T* later document published after the international filing date or priority date and not in conflict with the application but aimed to understand the principle or theory underlying the invention  |
| *A* document defining the general state of the art which is not considered to be of particular relevance  | *X* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone   |
| *E* earlier document published on or after the international filing date  | *Y* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art |
| *L* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified) | *A* document member of the same patent family  |
| *O* document referring to an oral disclosure, use, exhibition or other means  |  |
| *P* document published prior to the international filing date but later than the priority date claimed  |  |

Date of the actual completion of the international search

08 MAY 2000

Date of mailing of the international search report

12 JUN 2000

Name and mailing address of the ISA/US  
Commissioner of Patents and Trademarks  
Box PCT  
Washington, D.C. 20231

Facsimile No. (703) 305-3230

Authorized officer *James L. Grun*  
JAMES L. GRUN, PH.D.

Telephone No. (703) 308-0196

Form PCT/ISA/210 (second sheet) (July 1998)\*

## INTERNATIONAL SEARCH REPORT

International application No.  
PCT/US00/00855

## C (Continuation). DOCUMENTS CONSIDERED TO BE RELEVANT

| Category* | Citation of document, with indication, where appropriate, of the relevant passages  | Relevant to claim No. |
|-----------|---|-----------------------|
| Y         | GAO et al. Immunochemiluminometric Assay with Two Monoclonal Antibodies Against the N-Terminal Sequence of Human Parathyroid Hormone. Clinica Chimica Acta. 1996, Vol. 245, pages 39-59, see entire document. | 1-46                  |